

# EZMeta<sup>™</sup> Fatty Acid Synthetase (FAS) Colorimetric Activity Kit

Cat #: D-AKC2250 Size: 48T / 96T Storage: Stored at -20°C for 6 months

### **Product Information**

Detection range: 0.078-5 mmol/L (The detection range corresponds to the standard, and the actual content of sample is 0.156-10 mmol/L) Sensitivity: 0.039 mmol/L (The sensitivity corresponds to the standard, and the actual content of sample is 0.078 mmol/L) Applicable samples: Serum

Applicable samples: Serum

#### **Assay Principle**

High Density Lipoprotein Cholesterol (HDL-C) as one of the serum proteins, mainly synthesized by the liver. It's carrying cholesterol from surrounding tissues, which is then converted into bile acids or directly excreted from the intestine through bile. It is an anti-atherosclerotic plasma lipoprotein, a protective factor for coronary heart disease (CHD), commonly known as "vascular scavenger", and is an important reference index for the clinical diagnosis of coronary heart disease (CHD). EZMeta<sup>TM</sup> Fatty Acid Synthetase (FAS) Colorimetric Activity Kit provides a convenient tool for detection of HDL-C. The principle is that HDL-C in serum was separated with precipitant, and cholesterol esterase was used to catalyze the hydrolysis of cholesterol ester into free cholesterol (FC) and free fatty acid (FFA), so as to convert cholesterol ester into FC. Furthermore, cholesterol oxidase catalyzed FC oxidation to produce  $\Delta$ 4-cholestenone and  $H_2O_2$ . Then peroxidase catalyzes the oxidation of 4-amino-antipyrine and phenol with  $H_2O_2$  to produce red quinone compounds has a characteristic absorption peak at 500 nm.

**Biograde**tech Inc • 123 Beaver Street, Suite 204, Waltham, MA 02452 • USA E-mail : order@biogradetech.com • Phone: (800) 832-2612 • Website: www. biogradetech.com





### **Materials Supplied and Storage Conditions**

Kit components	Size		o
	48 T	96 T	Storage conditions
Extraction Buffer	10 mL	20 mL	4°C
Chromogen	7.5 mL	15 mL	4°C, protected from light
Enzyme Mix	1	1	-20°C, protected from light
Assay Buffer	5 mL	10 mL	4°C
Standard	0.25 mL	0.5 mL	-20°C, protected from light
Standard Diluent	5 mL	10 mL	4°C

#### **Materials Required but Not Supplied**

·Microplate reader or visible spectrophotometer capable of measuring absorbance at 500 nm

·96-well plate or microglass cuvette, precision pipettes, disposable pipette tips

·Refrigerated centrifuge, incubator, ice maker

·Deionized water

·Homogenizer (for tissue samples)

#### **Reagent Preparation**

Extraction Buffer: Ready to use as supplied. Equilibrate to room temperature before use. Store at 4°C.

**Chromogen:** Ready to use as supplied. Equilibrate to room temperature before use. Store at 4°C, protected from light.

Assay Buffer: Ready to use as supplied. Equilibrate to room temperature before use. Store at 4°C.

Standard Diluent: Ready to use as supplied. Equilibrate to room temperature before use. Store at 4°C.

Enzyme Mix: Add 5 mL Assay Buffer for 96 T or 2.5 mL Assay Buffer for 48 T to fully dissolve before use. This solution

can be stored at 4°C for one week or be stored at -20°C after aliquoting to avoid repeated freezing and thawing.

Standard: Containing 5 mmol/L cholesterol Standard. Store at -20°C, protected from light.

Working Reagent: Prepare 200  $\mu$ L Working Reagent for one well, add 50  $\mu$ L of dissolved Enzyme Mix and 150  $\mu$ L

Chromogen. Prepare Work Reagent before use and depend on your need.





Standard curve setting: Dilute 5 mmol/L cholesterol Standard with Standard Diluent to 5, 2.5, 1.25, 0.625, 0.313, 0.156,

0.078, 0 mmol/L Standard solution as shown in the table below.

Num.	Values of Standard	Volume of Standard Diluent	The concentration of Standard
	Volume of Standard	(μL)	(mmol/L)
Std.1	100 μL 5 mmol/L	0	5
Std.2	50 μL of Std.1 (5 mmol/L)	50	2.5
Std.3	50 μL of Std.2 (2.5 mmol/L)	50	1.25
Std.4	50 μL of Std.3 (1.25 mmol/L)	50	0.625
Std.5	50 µL of Std.4 (0.625	50	0.313
	mmol/L)	50	0.313
Std.6	50 µL of Std.5 (0.313	50	0.156
	mmol/L)	50	
Std.7	50 µL of Std.6 (0.156	50	0.078
	mmol/L)	00	0.078
Std.8	0	100	0

Note: Always prepare a fresh set of standards per use.

#### **Sample Preparation**

#### Note: Fresh samples are recommended. If not assayed immediately, samples can be stored at -80°C for one month.

Serum was separated within 3 h after blood collection and allow samples to clot for 30 min at room temperature before centrifugation for 15 min at 1,000 g. Remove serum and assay immediately or aliquot and store samples at -80°C. Avoid repeated freeze-thaw cycles. Thoroughly mix the serum sample with the Extraction Buffer in a ratio of 1:1 and let stand at 25°C for 15 min, then centrifuge at 2,000 g for 15 min at room temperature. Use supernatant for assay, and place it on ice to be tested.



### Assay Procedure

1. Preheat the microplate teader or visible spectrophotometer for more than 30 min, and adjust the wavelength to

500 nm, visible spectrophotometer was returned to zero with deionized water. Preheat the incubator to 37°C.

2. Add the following reagents respectively to the 96-well microplate or microglass cuvette:

Reagent	Standard Well (μL)	Test Well (μL)
Different Concentration Std.	20	0
Sample	0	20
Working Reagent	200	200

3. Mix well and incubate at 37°C for 30 min. Then reading the absorption values at 500 nm as A.

### **Data Analysis**

1. Drawing of standard curve

Each standard and sample absorption values subtract the zero standard (Std.8) absorption value to obtain  $\Delta A$ .

With the concentration of the standard solution as the y-axis and the  $\Delta A_{Standard}$  as the x-axis, draw the standard curve.

Substitute the  $\Delta A_{Test}$  into the equation to obtain the y value (mmol/L).

- 2. Calculate the content of HDL-C in sample
- HDL-C (mmol/L)=2×y
- Where: 2: Diluted multiples during sample preparation=(1+1)/1=2.

## **Typical Data**

Typical standard curve:





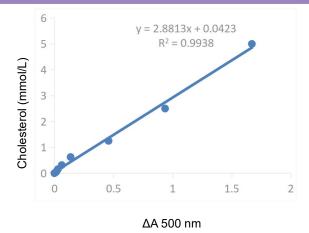


Figure 1. Standard curve for cholesterol.

### Disclaimer

The reagent is only used in the field of scientific research, not suitable for clinical diagnosis or other purposes.

